

FlexAssure™

FLEXIBLE INSTRUMENTS

Flexibility. Safety.
Strength. FlexAssure.




OmniGuide's FlexAssure flexible instruments with PolyBend technology have been successfully used in tens of thousands of intra-abdominal procedures without breakage because of its core features:

FLEXIBILITY: FlexAssure flexible instruments are the only CO₂ laser delivery instruments on the market constructed with a polymer core. Polymer plus multiple performance layers provide strength and durability while also increasing flexibility.

SAFETY: FlexAssure flexible instruments are the only CO₂ laser delivery instruments on the market engineered to safely contain laser energy in the event of a failure.

STRENGTH OF DESIGN: Stronger is superior. FlexAssure flexible instruments are the only CO₂ laser delivery instruments on the market that are fortified by a stainless steel jacket, ensuring that the instrument will not have shards break off into the patient. Other competitors build their CO₂ delivery instruments purely out of glass (silica) with no safety jacket design feature.

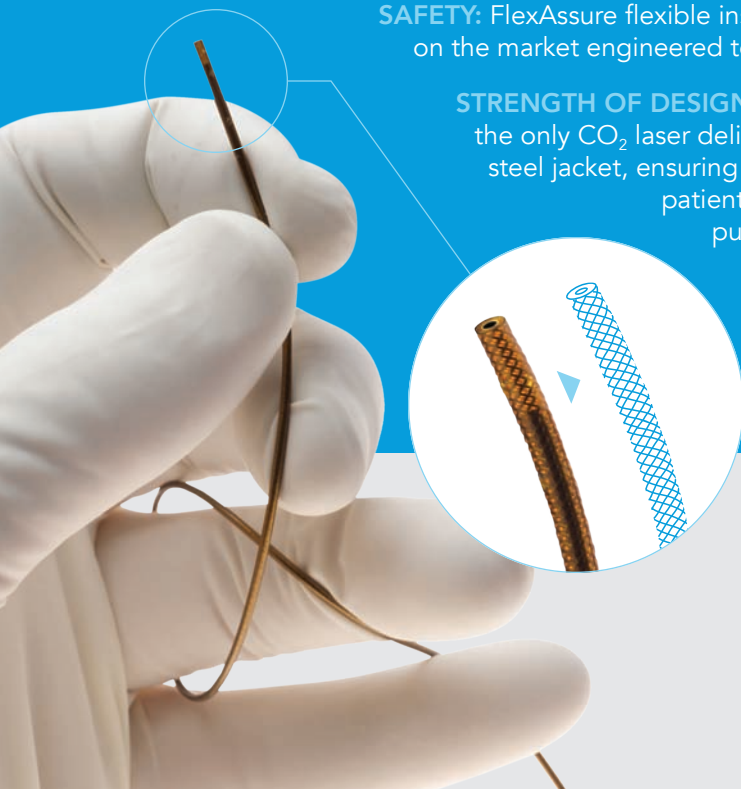

◀ Close up of a FlexAssure instrument highlighting the flexible stainless steel jacket.



OmniGuide's FlexGuide™ ULTRA instrument brings the precision of CO₂ energy to robotic surgery. FlexGuide ULTRA is indicated for use in robot-assisted procedures only, in conjunction with the BeamPath® ROBOTIC fibers to enable the surgeon to perform incision, excision, ablation, vaporization, and coagulation of body soft tissue.

Coiled Design

Flexibility to accommodate full robotic wrist articulation

Compatible with 5mm & 8mm Large Needle Drivers

No additional robotic arm or instrumentation required



Distal Fiber Stop


For easy loading and ensuring proper seating of fiber

4mm Spatula Tip

Efficient blunt dissection & convenient estimation of distance from the targeted tissue



With the ability to cut, ablate, and coagulate, the FlexGuide ULTRA and BeamPath Fibers can be used to address gynecologic pathologies such as endometriosis, myoma, adhesions and cysts.

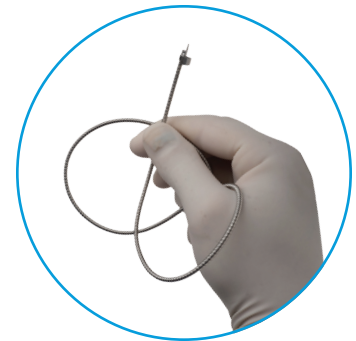

SKU	Description
420060-AA	BEACON Advanced Energy Laser System, 60 watt
420030-AA	BEACON Advanced Energy Laser System, 30 watt
BP-Robotic	BeamPath Robotic Fiber
GYN-HS-MIS	GYN MIS Robotic & Laparoscopic Inst Set 3HP + Tray
FLEX-ULTRA SET	FlexGuide ULTRA with Tray 1HP + Tray

For more information call 888-666-4484 or visit www.omni-guide.com

OmniGuide and OmniGuide Surgical are registered trademarks of OmniGuide, Inc. All content is copyright of OmniGuide, Inc. © 2019 See www.omni-guide.com/patents for information about the coverage of this OmniGuide product by U.S. and/or foreign patent rights. 1-0168-080-00-00 rev 0. To request copies use: MKT-024

888.666.4484 | 617.551.8444
customerservice@omni-guide.com
www.omni-guide.com

